



Non-ablative
resurfacing,
redefined.

CYNOSURE[®]
LUTRONIC[®]





Achieve efficacious results with our *Integrated Cooling System*

The MOSAIC integrated cooling system enables providers to safely and effectively address sensitive treatment areas around the eyes and mouth. Integrated air-cooling flows through the handpiece and onto the unique cold-compression sapphire tips, which cool and compress the skin. This synergy reduces pain, erythema, and edema, and is designed to enhance treatment efficacy, safety, and comfort.



Cold Compression Sapphire Tip



Handpiece Air-Cooling in Action







Adding new dimensions with *3D Controlled Chaos Technology™*

3D Controlled Chaos Technology (3D CCT) uses a machine learning model to add a new dimension to our patented CCT scanning.

By optimizing the distribution and depth of multi-depth microbeams, 3D CCT helps reduce the risk of bulk thermal heating, allowing for an increase in treatment coagulation volume per stamp, without sacrificing safety.

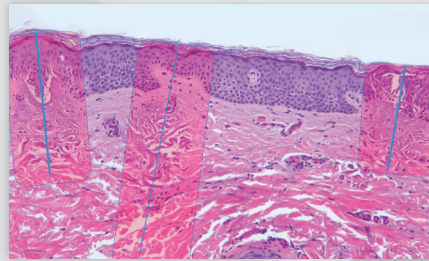

3D CCT Histology

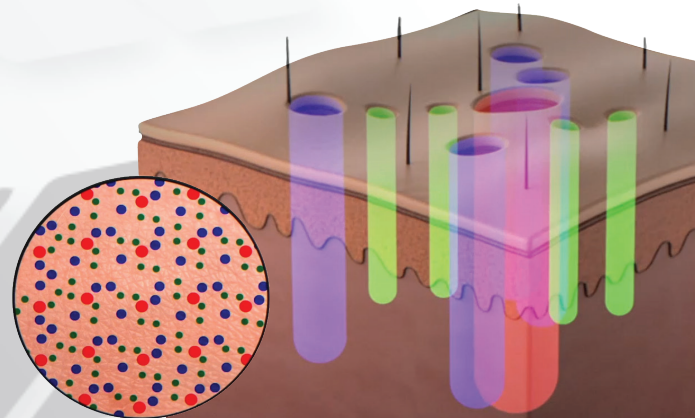


Real patients. Genuine concerns. *Life-changing results*✓

Courtesy of Robert Murgia, DO



Before and after 2 treatments

Courtesy of Amy Taub, MD

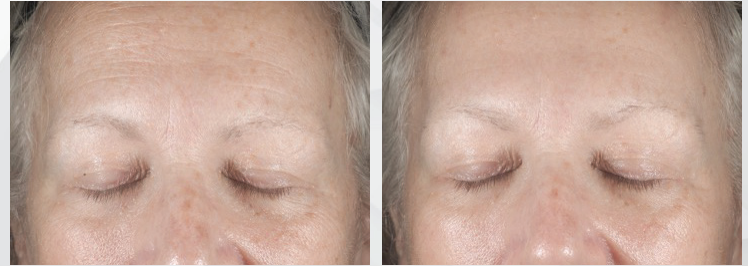

Before and after 2 treatments

Courtesy of Robert Murgia, DO



Before and after 2 treatments

Courtesy of Amy Taub, MD



Before and after 2 treatments



Distinguish your *services*

MOSAIC creates a new canvas in the skin by enhancing neocollagenesis and neoelastinogenesis in the mid-dermis, leading to revitalized and brightened skin. With its upgraded metal rollers and new sapphire tips, MOSAIC has unified rolling and stamping into one laser platform. Easy, air-cooled rolling enables efficient coverage and unmatched efficacy, while 3D CCT stamping remodels collagen in challenging focal areas.

MOSAIC is a highly controllable and customizable fractional laser platform; its flexible treatment range enables ever-changing utility in your ever-expanding practice.





Patient and practice *benefits*



3D CCT effectively stimulates mid-dermal coagulation, while sparing thermal injury



Reduced downtime with enhanced results



High quality textural improvement, especially around the eyes & mouth



Can be used in a complementary manner with other treatment modalities



Cold compression sapphire tips and air-cooling for a better patient experience



Comprehensive stand-alone resurfacing for high ROI





We're your trusted *partner on this journey/*

Technology

Forward thinking science
Clinically valuable innovations and capabilities

Product

Significant clinical evidence
Highest quality technological improvements

Practice Development

Industry expert consultants
Focused on growing your business

Clinical

Comprehensive virtual training
Personalized in-office training for beginner
and advanced teams

Service

Experienced technical support team 24/7
Expert field service team

Amplified Marketing Practice Support

Dedicated patient support team
Comprehensive turnkey marketing





Visit/
www.lutronic.com

Follow us:

@cynosurelasers

Cynosure, LLC

@lutronic_us

Lutronic

Intended for distribution in US only.

MOSAIC is indicated for the following:
dermatological procedures requiring the coagulation of soft tissue, as well as for skin
resurfacing procedures.

Cynosure, LLC and Lutronic Aesthetic, Inc. do not diagnose or treat medical conditions,
prescribe treatment protocols, provide medical advice, or otherwise practice medicine.
Patient should consult with a qualified healthcare professional for diagnosis and proper
treatment protocols/parameters to determine if treatment with any product is right for
them. Qualified healthcare professionals are responsible for diagnosis and treatment
protocols using products.

© 2024 CYNOSURE LUTRONIC, INC. All Rights Reserved. Mosaic, Controlled Chaos
Technology and 3D CCT are trademarks of Lutronic Aesthetics. PRD-1425
921-7073-001 R01



Attributes & Specifications

Attributes	Specifications
Laser Type	1550nm, Er: Glass
Laser Delivery System	Optical Fiber with Handpiece
Pulse Energy	4 to 70mJ
Pulse Rate	up to 447.4Hz
Cooling System	Integrated Air Cooling
Cold Compression Sapphire Tips	Circle: up to 12mm diameter Tx area Rectangle: up to 8x20mm Tx area Square: up to 20x20mm Tx area
Controlled Chaos Technology™	Exclusively available in 2D and 3D
Dimensions & Weight	16" W x 23" D x 69" H & 73lbs (33kg)
Electrical Requirements	AC100-240V, 50Hz/60Hz (Fuse AC125V/6.3A) Max Power: 350VA

